



PRINT DATE: Thursday, 25 August 2016



▶ VIEW DATASHEET ON-LINE www.chinesport.com/go/08390/

Status: Replaced

08390 FISIOTEK 2000 N

A direct-current motor, controlled by a microprocessor electronic board, operates the other Fisiotek 2000 units. All functions can be programmed by using the on-board control panel. START and STOP keys are also available on a hard-held control, so that patients can start and stop movement if needed.

The unit can be adjusted to fit the patient according to the limb length, as measured from the greater trochanter to the malleolus.

This device allows the ${\bf knee}$ and ${\bf hip}$ to be mobilized through flexion-extension movements.

It can accommodate limbs that are 72 - 100 cm long, which is equivalent to a patient's height of about 150-200 cm, respectively.

Programmable functions:

- knee flexion-extension range from 0° to 110°;
- speed from 0.8 to $4.6^{\circ}/$ sec;
- force 0-40 kg;
- session time;
- pause during flexion.

Product changes

Replaced

Warning! This article has been replaced with the article: XRI005





08390 - FISIOTEK 2000N





Technical specifications and dimensions

Mobilization :	Hip - Knee
Use of memory card :	-
Speed during flexo/extension control:	-
Speed control:	X
Force control :	X
Work time control:	X
Automatic extension increase :	-
Automatic flexion increase :	-
Pausa during extension :	-
Pausa during flexion :	X
Warm up cycles :	-
Knee movement breadth:	-5° - 110°
Hip movement breadth (mid limb):	10° - 70°
Speed:	0,8°/sec - 4,6°/sec.
Force:	0 - 40 kg
Power supply :	85 - 260 V; 50 -60 Hz
Electrical safety :	Class 1 B Standard EN 60601-
Electromagnetic compatibility :	Group 1 Class B Standard EN 60601-1-2
Classification according to EEC Directive 93/42:	Class II a
Weight:	14 kg

Documentation available online

DOCUMENTS IN ENGLISH

Brochure - Fisiotek 2000 Series .pdf 336.48 KB Insights - Fisiotek 2000 Series .pdf 8.24 KB User and Maintenance Manual - 08390 - User manual .pdf 1.4 MB

DOCUMENTS IN ITALIANO

Insights - Scheda tecnica .pdf 24.34 KB User and Maintenance Manual - 08390 - User manual .pdf 1.55 MB

DOCUMENTS IN ALL LANGUAGES

Declaration of conformity - 08390 - $FISIOTEK\ 2000\ N$ -▲ .pdf 213.89 KB Declaration of Conformity

⚠ The document is available for download in a reserved area

www.chinesport.com

FISIOTEK 2000 N



Accessories



01840 ACCESSORY FOR LIMBS SMALLER THAN 72 CM Rehabilitation exercises are possible for 61 - 72 cm long limbs with a 0° - 110° range .




01841 FISIOTEK TROLLEY Simple and functional, fi tted with non-slip supports for Fisiotek anda service shelf.


FISIOTEK 2000 N



More pictures of product

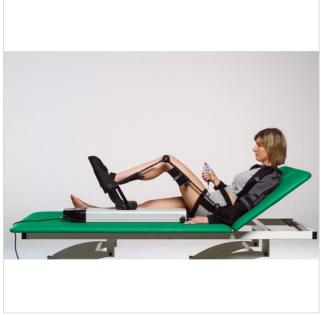

08390 - FISIOTEK 2000N

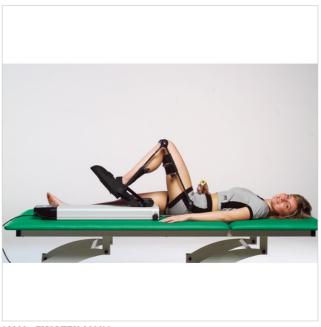

08390 - FISIOTEK 2000N